


# General Bed 35000EC1



### General Bed 35000EC1

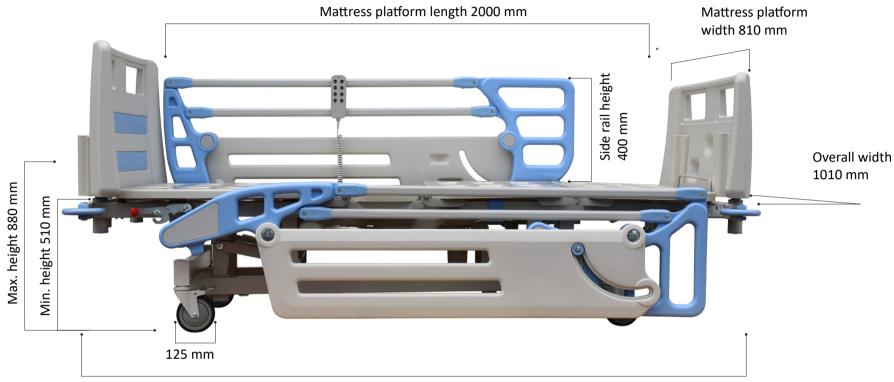

Overall length 2250 mm

### Technical Specifications

$2250\ \text{mm}\times1010\ \text{mm}$
$2000 \text{ mm} \times 810 \text{ mm}$
880 mm - 510 mm
400 mm
. 0 - 80°
0 - 30°
135 kg
250 kg


Casters	. 125 mm
Power input max.	250 VA
Main power system	220V – 50 Hz, 24 V DC
Electrical current protection	Class I
Electrical shock protection	Туре В
Operation mode	10%, max.: 2 min/18 min
IP rating	IP 66

Technical and Quality Assurance Standards IEC 60601 - 2 - 52 ISO 13485 ISO 14971

## Positions





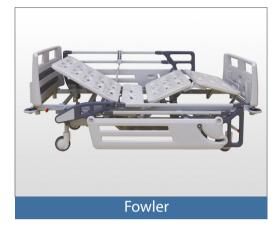







### General Bed 35000EC1



#### Bed Ends

The high quality, seamless blow molded bed ends of model 35000EC1 are lightweight, durable, and easily washable and steam cleaned. In case of emergency, the caregiver can easily remove the bed ends by grasping the hand holds and lifting up, offering better access to the patient's head and upper body. Additionally, the bed end's integrated steel fittings further enforce the sturdiness of the structure and help to ensure smooth, reliable transportation.



#### Pendant Control

The user friendly handset control offers an array of illustrated push buttons, making adjusting the position of the bed simple. The pendant controller can be placed on both of the side rails for easy access.



#### Side Rails


The parallel side rails of the 35000EC1 offer a unique and stylish design produced from a combination of blow molded plastic and steel bars to provide a protective height and length. The side rails adhere to the highest standards of safety to prevent patient falls and entrapment."

### Features



#### Mattress Platforms

The ABS mattress platforms are produced by injection molding. The platforms boast a number of helpful features, like smooth surfaces, high durability, and a lightweight, easily transportable frame enabling quick and painless removal for cleaning and disinfection, including steam cleaning.

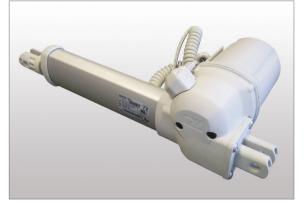

#### Actuators

The powerful actuators of the 35000EC1 maneuver the bed seamlessly and quietly into a variety of positions. With an IP66 rating, the actuators are dust and water resistant and are easily accessible for cleaning and disinfection.

#### Bumpers

Rotary on all f to miti during

Rotary bumpers are located on all four corners of the bed to mitigate potential damage during unintended collisions



#### Casters

The 130 mm casters move smoothly, ensuring quick and easy transport for the nursing staff. The diagonal locks can be activated to prevent unintended movements.



Urine Bag Holder Urine bag holders can be easily installed in one step on either side of the bed.